


Icon I3 Triple Layer Medium Risk Pressure Relieving Mattress

The Icon I3 support surface range provides true leadership in static pressure care in a preventative setting. Tried and trusted in facilities nation-wide, the Icon I3 range combines outstanding durability, a high standard and consistency of care and excellent value. An extensive variety of configurations are available providing for a wide range of patient needs and requirements.

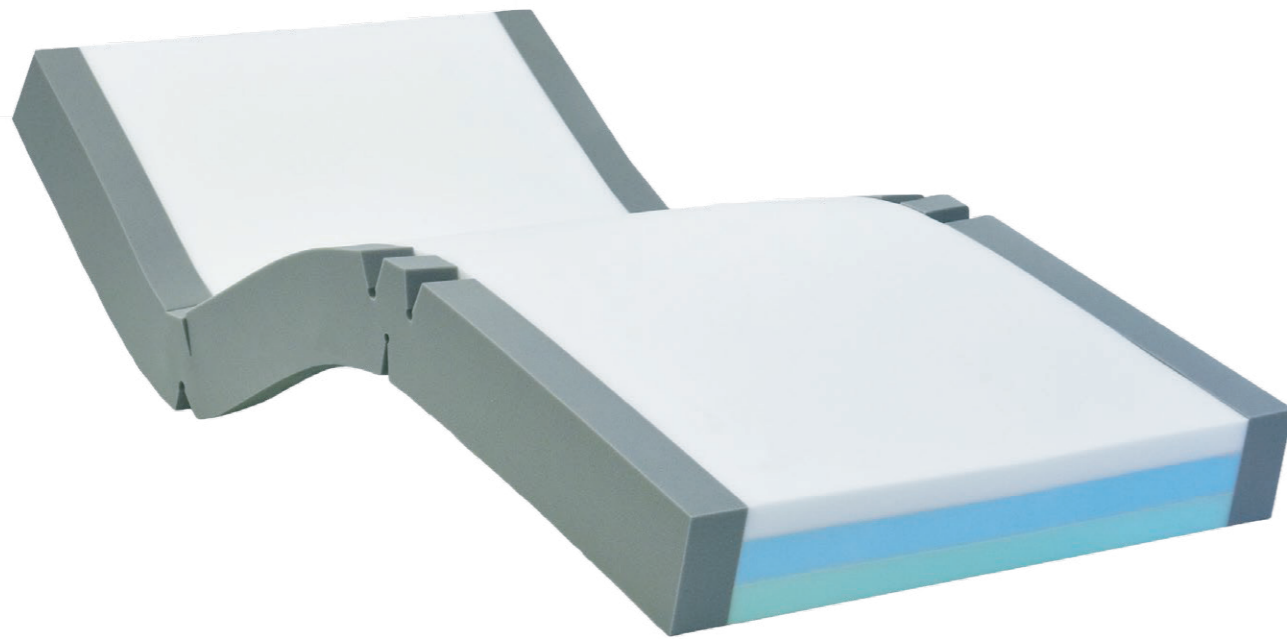

/ Icon I3


CODE	SIZE (MM)	
I3S1000	1980x880x150	Single
I3SL1000	2030x880x150	Single Long
I3KS1000	1980x1050x150	King Single
I3KSL1000	2030x1050x150	King Single Long
I3EW1000	1980x1150x150	Extra-Wide
I3EWL1000	2030x1150x150	Extra Wide Long
I3D1000	1980x1350x150	Double
I3DL1000	2030x1350x150	Double Long
I3Q1000	1980x1530x150	Queen
I3QL1000	2030x1530x150	Queen Long
THERAPEUTIC SUPPORT LOADING	40-150KG	
APPLICATIONS	✓ HOSPITAL WARDS ✓ AGED CARE ✓ PALLIATIVE CARE ✓ HOME CARE	
WARRANTY	5-YEAR FOAM CORE WARRANTY 2-YEAR COVER WARRANTY	

Features




AIROFORM MEMORY IMMERSION INTERFACE LAYER

- ✓ Alleviates localised pressure points
- ✓ Provides maximum contact with patient for distribution of interface pressure
- ✓ High temperature stability to deliver consistent immersion, performance and comfort
- ✓ Visco-elastic cell structure provides low-shear surface
- ✓ Conforms to EPUAP, NPUAP & PPIA Guidelines for quality and comfort

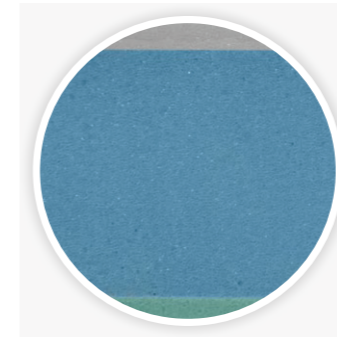


DELUXE GRADE MID-LAYER

- ✓ Open cell high quality foam
- ✓ Provides exceptional comfort and postural support for repositioning and patient care

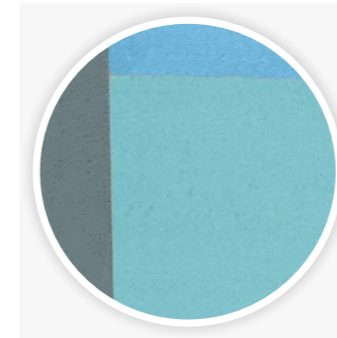

FIRM FOUNDATION LAYER

- ✓ Provides protection against 'bottoming-out' during transfer, repositioning and whilst in the inclined/fowler position
- ✓ Provides exceptional comfort, support for repositioning and patient care

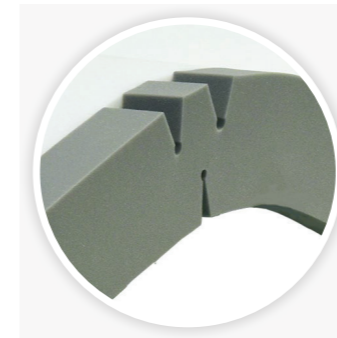

PROFILE STABILISING SIDES

- ✓ Enables ease of primary patient care
- ✓ Facilitates patient transfers
- ✓ Unique hinge points allow mattress to easily conform to varied bed positions
- ✓ Encourages central patient positioning

Covers



NEW! LUXURY QUILTED COTTON COVER

- ✓ Available on all Forté Mattresses
- ✓ Contact Forté for further information

/ Mattress & Cushion Covers

	ADVANTIFLEX CLASSIC	PREMIFLEX CLASSIC	PREMIFLEX ULTRA
High Frequency Welded Seams	✓	✓	✓
Integrated Zipper	✓	✓	✓
Waterfall Zipper Protection Flap	✓	✓	✓
4-Way Stretch Patient Contact Surface	✓	✓	✓
Cover Warranty	2 Years	4 Years	4 Years
Chlorine / Sodium Hypochlorite Resistance	1,000PPM	10,000PPM	10,000PPM
Heavy Duty PVC1100 Base	✓		
Heavy Duty PVC1200 Base		✓	
Heavy Duty DoubleCoated Polyurethane Base			✓
Premium Comfort, Non-Abrasive and Durable Upper	✓		
Superior Durability, Increased Shear Protection, Ultimate Comfort Layer		✓	✓
Full Vapour Permeability	✓	✓	✓

Features

Advantiflex Classic Pressure Care Cover

Belgium made polyurethane coated upper and heavy duty PVC Base ensures that Forté Advantiflex Covers deliver consistent care, comfort and durability. Utilising premium quality medical grade fabrics, Forté Combination Covers are a tried and proven medical grade Pressure Relieving Cover.

Premiflex Classic Pressure Care Cover

Forté Premiflex Classic Covers provide outstanding infection control, hygiene and high performance within the Acute Care Setting. Manufactured with a High-Quality Belgium Halogen Free Polyurethane 4-Way Stretch, Moisture Vapour Permeable Anti-Shearing Interface Material, the Premiflex Cover is an effective Pressure Injury risk mitigation mattress cover.


Premiflex Ultra Pressure Care Cover

Forté Premiflex Ultra Covers provide the ultimate protection and performance for mattresses in an acute care setting. Superior pressure injury protection is achieved through 4-way stretch non-hammocking and anti-shearing characteristics, Premiflex Ultra Covers have a durable double coated polyurethane base that is moisture vapour permeable allowing foam longevity by allowing the core to breath.




HIGH FREQUENCY SEALED AND WELDED SEAMS

- ✓ Electromagnetically fused seams protecting against risk of penetration of fluid
- ✓ Meets and exceeds infection control requirements
- ✓ Reinforced with additional heavy duty stitching at stress points
- ✓ 100% sealed to prevent ingress of fluid or infection

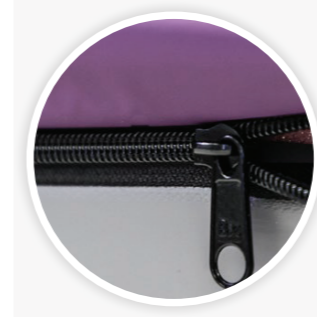

INTEGRATED ZIPPER

- ✓ Allows easy cover removal
- ✓ Robust and durable
- ✓ High frequency welded attachment for optimal protection against fluid ingress

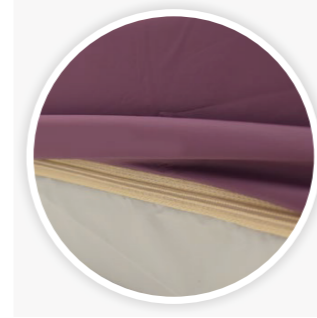

WATERFALL ZIPPER PROTECTION FLAP

- ✓ Protects against fluid ingress at point of zipper
- ✓ Meets infection control guidelines