
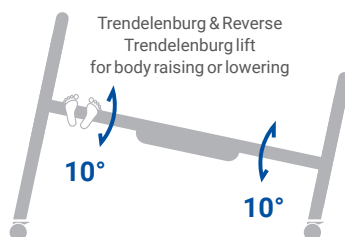
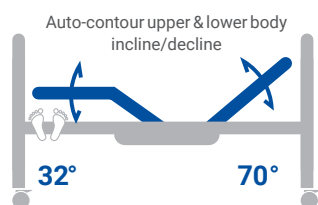
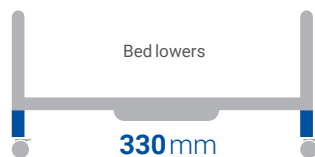
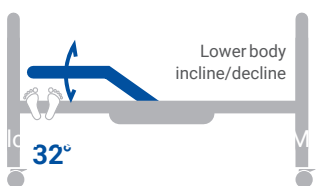
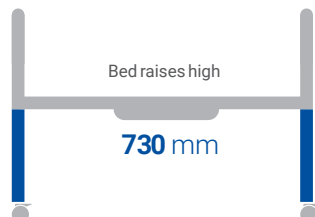
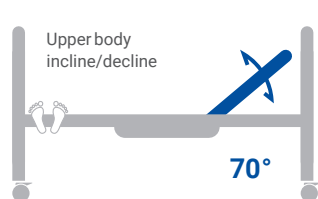


# BetterLiving™

## Community Bed

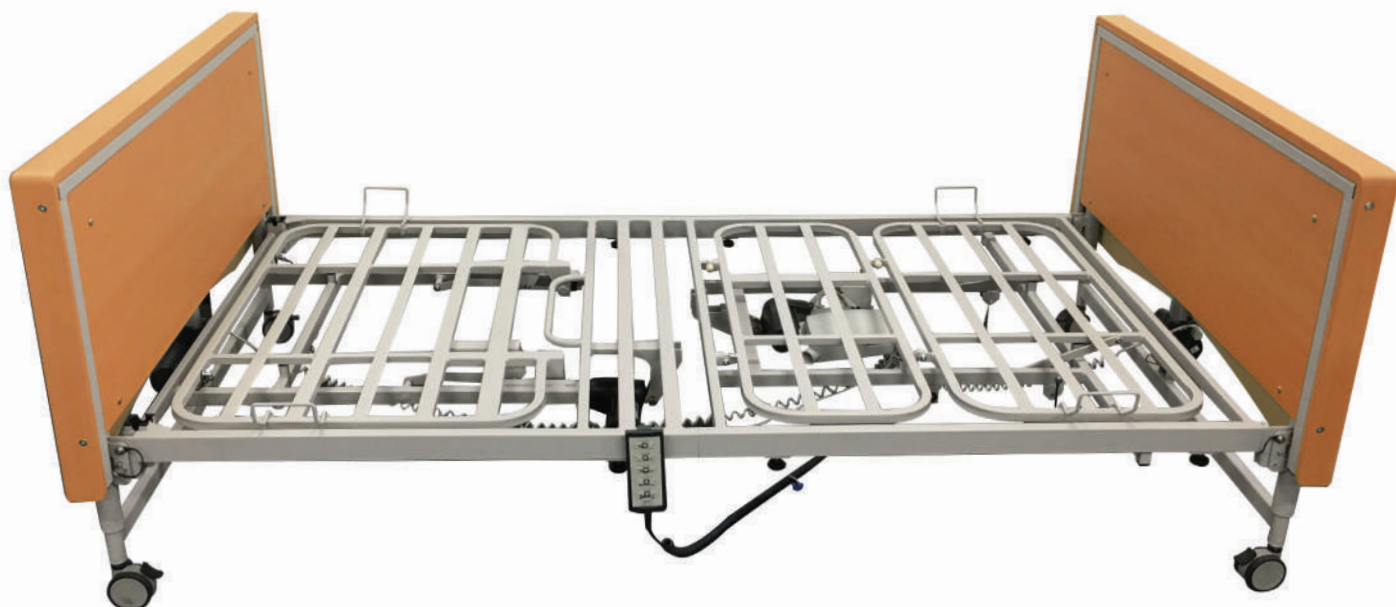

### Multiple Positioning and Reduced Shear

The BetterLiving Community Bed is incredibly practical yet versatile to suit a range of home care needs. Operates as a four part, profiling bed with electric high-low height adjustment, back rest, knee break and Trendelenburg tilt – to deliver multiple positioning options to suit every need.

The Auto Contour feature simultaneously raises both the head and knees, which combined with auto-regressive backrest helps prevent painful skin shearing.

### Quality Design

Powered by a high quality Linak electric drive and constructed from durable powder coated steel. The frame includes a 200mm Length Extension and removable mattress retainers.

Designed with built-in storage features for the Length Extension and Transport Bracket which minimises misplacement when not in use. Fitted with four high quality individually lockable castors for easy maneuverability.

### Easy, Safe Operation

Lightweight, water resistant hand controller with soft touch buttons for raising, lowering and positioning, including one touch auto contouring.

Convenient inbuilt hanging hook for easy storage and access. Features lockout function, ideal for patients incapable of operating the bed due to cognitive impairment.

Each function can be individually locked out using the key, which is required to unlock and reactivate functionality.



Phone 1300 734 862  
support@patienthandling.com.au

## Community Bed

SIZE	Single	King Single
CODE	BE1005	BE2005
SAFE WORKING LOAD	220 kg	220 kg
MAXIMUM PATIENT WEIGHT	180 kg	180 kg
WARRANTY	5 year bed frame / 2 year electrical / 1 year castor	5 year bed frame / 2 year electrical / 1 year castor
FRAME MATERIAL	Powder Coated Steel	Powder Coated Steel
DIMENSIONS (MM)	LENGTH	2170 to 2370 (Extended)
	WIDTH	995
HEIGHT RANGE	330 to 730 mm	330 to 730 mm
MATTRESS SIZE(MM)	LENGTH	1960
	WIDTH	865
RECOMMENDED	HEIGHT	170 maximum
BED FRAME WEIGHT (KG)*	74 kg	84 kg
BACKREST ANGLE	0 to 70°	0 to 70°
KNEE BREAK ANGLE	0 to 32° (4 Section Bed)	0 to 32° (4 Section Bed)
TRENDELENBURG / REVERSE TILT	0 to 10°	0 to 10°
COMPLIANCES	IEC 60601-1: 2005+A1: 2012 [Medical Electrical Equipment Part 1]	
	IEC 60601-2-52: 2012 [Medical Beds Safety and Essential Performance]	
	IEC 61000-3: 2008 [EMC Electromagnetic compatibility]	
	IEC 61000-4: 2009 [EMC Electromagnetic compatibility]	
ACCESSORIES	CISPR 11: 2004+A2: 2006, Group 1 Class B [EMC Electromagnetic compatibility]	



BE6003  
Side Rail, Three Quarter  
Length (Sold as Pair)



BE8003  
IV Pole



BE7003  
Self Help Pole



BE4017  
Transport Bracket  
(Sold as Pair)



BE9002  
Grab Bar



\*Shown in storage mode  
with optional transport  
bracket BE4017