

# Raizer II

## Lifting Chair

Raizer II is a mobile lifting chair that helps people who have fallen on the floor to get up again in just a few minutes. Raizer can be operated by a single assistant and requires only limited physical effort beyond the help of a supportive hand. Raizer II is an updated version - new design and features - of the original Raizer.

### Technical Specifications

MODEL	Raizer II
LIFTING CAPACITY	max. 150kg
LIFTING TIME	20-30 secs
WEIGHT	Seat: 8.5 kg / Legs and back supports: 4.5 kg / Total: 13 kg
NUMBER OF LIFTS WHEN FULLY CHARGED	40 at max. load. Approx. 80 with an average load
CHARGING TIME	Up to 6 hours
CHARGING OF EMPTY BATTERY	10-15 min. = 1 lift
CHARGERS	USB charger and 15W charger for standard power socket



PRODUCT CODE	PRODUCT NAME
MHLU-105029	RAIZER II Lifting Chair - including remote control, carry case, cover for seat, battery charger 15W, USB charger and packing.
MHLU-107464	Head Rest
MHLU-107273	Raizer II Trolley
MHLU-107250	Seatbelt
MHLU-103314	Triple Wheel for Trolley



# Raizer II

## Lifting Chair




### EVEN MORE CLEANING FRIENDLY

- All surfaces of Raizer II are designed to be very easy to clean
- New surface makes it possible to disinfect Raizer II




### EASIER TO ASSEMBLE

- 4 identical legs and 2 identical back rests ensure correct assembly every time
- Intuitive sound and LED indicators for correct assembly
- Improved carry bag and trolley

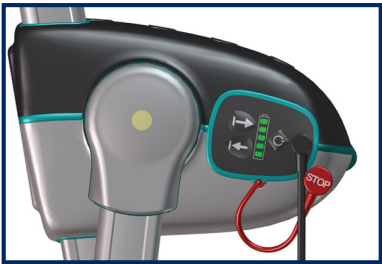

### NEW SEATBELT POSITION

- Seatbelt can be removed for washing
- Independently adjustable seat belt placement, on the back rest, ensures improved upper body stability when lifting fallen person

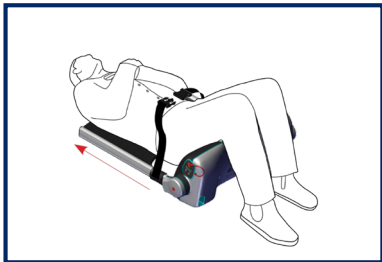

### REMOTE CONTROL AND EMERGENCY STOP

- Purpose built storage cavity for remote, moulded to the side of seat casing
- Emergency stop button is more intuitive
- Service reminder will indicate when Raizer II needs service
- Battery indicator



### END STOP ENSURES CORRECT POSITION

- End stop ensures correct positioning for use and disassembly